



EMPOWERMENT RULES.

With its new clean design and R-Net electronics, there is no protesting the fact the C500 VS is the best seat in town. Perfect for both indoor and outdoor use, the C500 VS has enough power and torque to navigate through rough terrain. Yet the independent suspension of the C500 VS provides a smooth ride. Users can maneuver in tight spaces with a 715 mm (28") turning radius, perfect for making your way through a crowd. Plus, the seating system with three programmable memory positions allows you to stand upright from any position as well as to drive while standing.

C500
VS JR/SR STANDER

POWER FUNCTIONS



R-Net (Standard).
Alternative Drive Controls
available. (Option)



Power sit to stand from any position. Drive at reduced speed in standing position.

SHROUD COLORS



Arctic White



Crystal Blue



Vintage Wine



Metallic Green



Platinum Silver



Black Onyx




Power adjustable recline
operates independently from
power adjustable legrests.



Power adjustable tilt to 25°.



Power adjustable legrests
operate independently from
power adjustable recline.

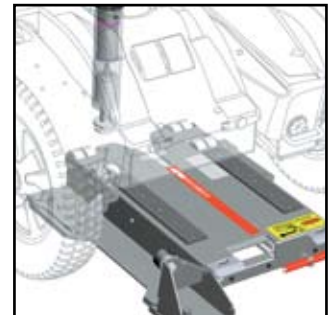

Permobil C: Lock your Permobil
into position for safer travel by car.
(Option)

C500 VS™ JR/SR SPECIFICATIONS

Weight incl. Batteries	171 kg (377 lbs)	Max Speed Standard	7.5 km/h (4.7 mph)
Base length	1215 mm (47.8")	Max Speed Option S model	10 km/h (6.5 mph)
Base width	675 mm (26.6")	Driving speed standing	2 km/h (1.3 mph)
Turning radius	715 mm (28")	Driving speed tilted	½ speed at 135° backrest angle
Seat to floor height	470–535 mm (18.5"–21")	Range	30-40 km (18-25 miles)
Seat elevator (option)	250 mm max. (10")	Battery Capacity	2 x 73 Ah
Power tilt angle (option)	25° max.	Max user weight	Jr 70 kg (155 lbs) Sr 120 kg (265 lbs)
Power recline (option)	180° max.	Max user weight stand and drive	Jr 70 kg (155 lbs) Sr 100 kg (220 lbs)
Seat Widths	430/480 mm (17"/19")		
Seat Depths	Jr 360 mm (14") Sr 420–550 mm (16.5"–21.5")		
Backrest heights	550/620 mm (21.5"/24.5")		



THE POWER OF MOBILITY

Permobil AB.
Box 120, SE 861 23 Timrå, Sweden.
Phone: +46 60 59 59 00.
www.permobil.com