



BRILLIANTLY SIMPLE


Now we have developed the perfect combination of daring design and robust strength. We have created a new generation of seats for our two compact frames C300 and C350. The C300TS is a front-wheel drive Permobil that offers everything you expect of a reliable wheelchair, yet is priced to suit the pocketbook. The simple concept of the TS seat gives it comfortable form. The seat frame has integrated UniTrack rails that you may use to adjust seat depth and arm rest height. The rails also may be used to position any attached supports.

C300
TS

POWER FUNCTIONS



FRAME COLOUR

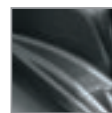

Black

Powered adjustment of back rest and seat inclination



VR2 control with four optional functions



Continuously adjustable seat depth and arm rest height



Divided rims in same colour as frame



The back rest and seat are made of closely-spaced, extra strong Velcro tape, so the seat can be adjusted to suit the body shape



Powered seat height adjustment



Powered leg rest adjustment

ISO 7176-19 / ISO 10542-3
EN12182/EN12184



SPECIFICATIONS

Length	108 cm	Legrest angle - adjustable	90° - 150°
Width	62 cm	Battery capacity	2 x 60 AH
Min. transport length	88,5 cm	Range	25 - 35 km
Seat height	51 cm	Max. speed	7 km/hr
Seat height with electric seat lift	51 - 71 cm	Turning radius (ISO 7176-5)	62,5 cm
Weight incl. batteries	128 kg	Obstacle capability	7 cm
Max. user weight	136 kg	Suspension	Yes
Seat width	38 - 43 - 48 - 53 cm	Electronics	VR-2
Seat depth	38 - 55 cm	Electric seat lift	Yes
Backrest width	38 - 43 - 48 - 53 cm	Electric tilt adjustment	Yes
Backrest height	45 - 60 cm	Electric legrest adjustment	Yes
Distance between armrests	38 - 53 cm	Electric backrest adjustment	Yes
Backrest angle - adjustable	90° - 125°	Warranty	1 Year



The Power of Mobility

Permobil Ltd
Oceanic, Ground Floor
Navigation Business Park
Waters Meeting Road
Bolton, BL1 8SW
info@permobil.co.uk

www.permobil.com