



standard
ESP Enhanced Steering Performance

MAKING THE GRADE

The C400 VS redefines the way standing is incorporated into a child's daily routine and allows the user to drive while in a standing position. Standing can be achieved from either a sitting or a horizontal position, enabling the child to move about throughout the day and remain comfortable. Designed specifically for children, the C400 VS features a front-wheel drive base and a rugged suspension system. The small footprint of the C400 VS combines maneuverability with advanced obstacle climbing ability. On top of which, the modular design of the Corpus seat ensures optimal fit with room for growth. With the C400 VS, every child can be top of the class.

C400 VS




POWER FUNCTIONS




R-Net electronics with ICS

Power sit to stand from any position.
Drive at reduced speed in standing position.

FRAME COLOURS

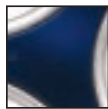

Dark blue metallic

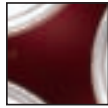

Red metallic



Green metallic



Black metallic

ISO 7176-19 / ISO 10542-3
EN12182/EN12184



Power adjustable recline



Power adjustable tilt



Power adjustable legrests operate independently from power adjustable recline

SPECIFICATIONS

Length	118,5 cm	Legrest angle - adjustable	90° - 180°
Width	62 cm	Battery capacity	2 x 73 AH
Min. transport length	105,5 cm	Range	30 - 35 km
Seat height	53 cm	Max. speed	7 km/hr
Seat height with electric seat lift	53 - 78 cm	Turning radius (ISO 7176-5)	66,5 cm
Weight incl. batteries	146 kg	Obstacle capability	6 cm
Max. user weight	70 - 100 - 120 kg	Suspension	Yes
Seat width	43 - 48 cm	Electronics	R-Net
Seat depth	41 - 60 cm	Electric seat lift	Yes
Backrest width	34 - 39 - 44 cm	Electric tilt adjustment	Yes
Backrest height	55 - 62 cm	Electric legrest adjustment	Yes
Distance between armrests	35 - 40 - 45 cm	Electric backrest adjustment	Yes
Backrest angle - adjustable	90° - 180°	Warranty	1 Year



The Power of Mobility

Permobil Ltd.
Oceanic, Ground Floor
Navigation Business Park
Waters Meeting Road
Bolton, BL1 8SW
info@permobil.co.uk

www.permobil.com