



KOALA MINIFLEX



MAXIMUM SPEED

6,5 KM/H



RANGE

20-25 KM



MOTOR POWER

2 X 300 WATTS



Freewheel lever



Easily accessible charging point



Pivoting rear wheels for responsive steering

The Power of Mobility

PAEDIATRIC WHEELCHAIRS

KOALA MINIFLEX

The Koala is the ideal wheelchair for children from toddlers to young teen-agers. It is fitted with the same robust motors and front wheel drive system as Permobil wheelchairs for adults but adapted to the needs of children.

The Miniflex seating system is a flexible seating system for growing children with special seating needs. There is a choice of three sizes and different body supports. Possible options include the fitting of a seating orthosis. When combined with an adult chassis, the Miniflex seating system can be easily exchanged for an adult seating system once the child has grown out of the Miniflex seating system.



Knee supports



Lateral supports



Thigh supports



The simple joystick controls are ideal for small children

TECHNICAL SPECIFICATIONS

LENGTH, MM	965	LEGREST ANGLE – ADJUSTABLE	90°-150°
WIDTH, MM	560	BATTERY CAPACITY , AH	2x40
MIN. TRANSPORT LENGTH, MM	710	CHARGING TIME, HOURS	8
SEAT HEIGHT, MM	410	RANGE, KM	20-25
SEAT HEIGHT WITH ELECTRIC SEAT LIFT, MM	410-660	MAX. SPEED, KM /H	6,5
WEIGHT INCL . BATERIES, KG	85	BRAKING DISTANCE , MM	1000
MAX . USER WEIGHT, KG	35	TURNING RADIUS PER ISO 7176-5	570
SEAT WIDTH, MM	250-300/350	OBSTACLE CAPABILITY , MM	60
SEAT DEPTH, MM	270/320/370	SUSPENSION	No
BACKREST HEIGHT, MM	325/375/425	ELECTRIC SEAT LIFT	Yes
ARMREST HEIGHT, MM	160-290	ELECTRIC TILT ADJUSTMENT	Yes
DISTANCE BETWEEN ARMRESTS, MM	215-375	ELECTRIC LEGREST ADJUSTMENT	No
DISTANCE BETWEEN SEAT CUSHION AND FOOTPLATE , MM	175-340	ELECTRIC BACKREST ADJUSTMENT	No
BACKREST ANGLE – ADJUSTABLE	90°-165°	CONTROL SYSTEMS	Yes
		WARRANTY	1 Year

ELECTRIC ADJUSTMENTS




ISO 7176-19 / ISO 10542-3
EN12182/EN12184



COLORES



YELLOW

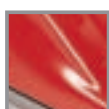

RED

Permobil Ltd.
Oceanic, Ground Floor
Navigation Business Park
Waters Meeting Road
Bolton, BL1 8SW
info@permobil.co.uk

www.permobil.com