



C500 RS



MAXIMUM SPEED

5 MPH




DRIVING RANGE

19-25 MILES



GROUND CLEARANCE

3 INCHES

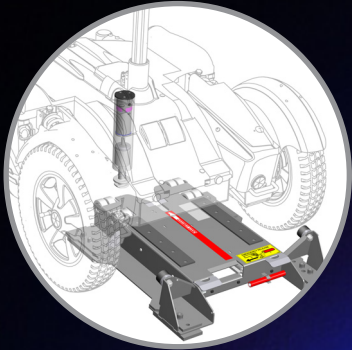

Permolock
Compatible



75° Total Motion Tilt
60° Posterior Tilt and
15° of Anterior Tilt



Dual Action Power
Recline



ICS - Intelligent Control
System. Program up to 3
memory seat functions; as
easy as tuning a car radio.



THE POWER OF MOBILITY

C500 RS



The all new C500 is built from timeless, classic Permobil values: superior quality, sophisticated technology and the greatest forms of functionality. Down to the smallest details, this base is exhaustingly refined and redesigned to be the best seat in its category. Its effortless performance tackles obstacles with unmatched power and increased maneuverability.

The RS (Rehab System) seat is "Inspired by nature and created by Permobil". The RS is the ultimate in modular seating systems. Featuring an integrated support system and Intelligent Control System (ICS), the RS offers a multitude of power seat functions that accommodate most, if not all, functional postures. The ICS allows for the utmost in customization with users being able to program up to 3 memory seat positions through the switch box. One push of a button will take you to your frequently used positions. Plus, the ICS offers synchronized movement of the actuators; VariSmart actuators offer real-time positioning and eliminate the need for mechanical limit switches. All this to allow users to achieve the most effective seating position desired...truly a divine creation.




Convenient Cup Holder and Cell Phone Holder



Adjustable Seat Width and Depth - 16" - 22" adjustment in 2" Increments



Adjustable Back Width - 14" - 22" Infinite Adjustment

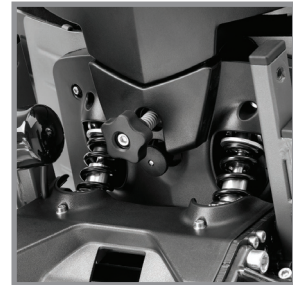

Back Rest Suspension - Adjustable twin shock absorbers provide a smooth comfortable ride.

TECHNICAL SPECIFICATIONS

LENGTH	44"	RANGE	19-25 Miles
WIDTH	25.75"-28.5"	MAX. SPEED, KM/H	5 MPH
MIN. TRANSPORT LENGTH	33"	TURNING RADIUS PER ISO 7176-5	43.75" Radius
SEAT TO FLOOR HEIGHT	18.5"*	GROUND CLEARANCE	3 Inches
POWER SEAT ELEVATOR	10"	SUSPENSION	Yes
WEIGHT INCL. BATTERIES	375 LBS	ALTERNATE DRIVE SYSTEMS	Yes,
MAX. USER WEIGHT	300 LBS	WARRANTY	2YR
SEAT WIDTH	16" - 22"		
SEAT DEPTH	16"-22"		
BACKREST HEIGHT	18.5", 22.5", 26.5"		
ARMREST HEIGHT	7.5"-14.5"		
DISTANCE BETWEEN ARMRESTS	16"-22" Infinite Adjustment		
ADJUSTABLE LEG LENGTH,	15.5" -22"		
BACKREST ANGLE - ADJUSTABLE	80°-165°		
LEGREST ANGLE - ADJUSTABLE	80°-150°		

POWER OPTIONS



COLOR OPTIONS



GREEN METALLIC

VINTAGE WINE

CRYSTAL BLUE

METALLIC CHARCOAL GRAY

SILVER/BLACK METALLIC

Specifications are subject to change

Permobil Inc.
6961 Eastgate Blvd.
Lebanon, TN 37090
Tel (800) 736-0925
Fax (800) 231-3256
info@permobilus.com
www.permobil.com